


## Active Transfer Device for Use in All Settings

Easy to use and non-powered, Hoyer Up is a manual stand-aid and seated transfer device suitable for sling-free transfers with caregiver assistance.

Hoyer Up is designed for users who have some weight-bearing ability and can participate in the transfer promoting patient rehabilitation and functional independence. The Up simplifies care in all settings and helps maintain more activities of daily living.



### UNIQUE MODULAR DESIGN

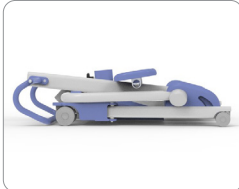


Simplifies inventory and transport with an efficient flat-pack design that requires no tools for assembly

### SWING SEAT




Secure padded swing seat surface delivers comfort and stability. Seats fold flat against the side arms when not in use.

### COMFORTABLE KNEE PAD

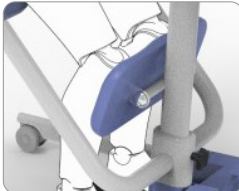

Enhances comfort and accommodates a larger range of patient sizes with a contoured, height adjustable knee pad.

### ADJUSTABLE BASE

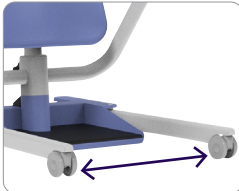

Increases patient care accessibility with a simple and robust mechanical foot pedal providing open/close functionality.



**SAFE WORKING LOAD: 440 LB / 200 KG**

### SPECIFICATION

Overall Length	37.4"	950 mm
Overall width	24.8"	630 mm
Overall Height	43.3"	1100 mm
Legs Open - External Width	33.3"	845 mm
Legs Open - Internal Width	28.5"	725 mm
Legs Closed - External Width	23.6"	600 mm
Legs Closed - Internal Width	18.9"	480 mm
Leg Height	4.3"	110 mm
Casters	3" Front	4" Rear

### WEIGHTS

Total	72.8 lb	33 kg
-------	---------	-------

Specifications may vary due to manufacturing tolerances.