


# ActiveCARE™ Mattress Series

---



**8 Interconnected Air Cells**

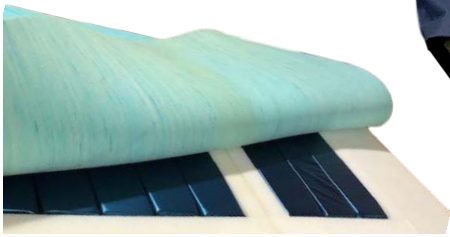

**Gel Infused Memory Foam**

## StatMatt Reactive Mattress


- **Zoned Reactive Air System** — eight horizontal foam filled air cells react to an individuals movement and weight. Interconnected bladders with pressure valves deflate, transfer air, and offer superior automatic pressure redistribution therapy intended to reduce risk of pressure injuries while improving comfort.
- **Visco Gel Memory Foam** - 2" gel infused memory foam topper contours to individuals body and increases comfort while improving envelopment for less interface pressure.
- **4-Way Stretch Cover**— when combined with memory foam and foam filled air bladders it provides superior immersion while reducing shear and friction plus and as an added benefit it has been treated with antimicrobial properties.
- **ActiveCare Bed Design** — multiple articulation cuts allow the mattress with air bladders to contour to the advanced positions the ActiveCare Bed supports such as full chair, high fowler position (almost 90 degrees) as well as additional bolsters to support egress points.
- **Separate Heel Pillow** — gentle heel slope with super soft foam heel pillow to offer superior pressure redistribution properties in the heel section.

## MATTRESS SPECIFICATIONS

**Mattress Weight Capacity**.....36" - 500lbs  
42" - 600lbs

**Mattress Specs**..... 36"(W) x 80" (L) x 7" (H)  
42"(W) x 80" (L) x 7" (H)

**Mattress Warranty**..... 1 year

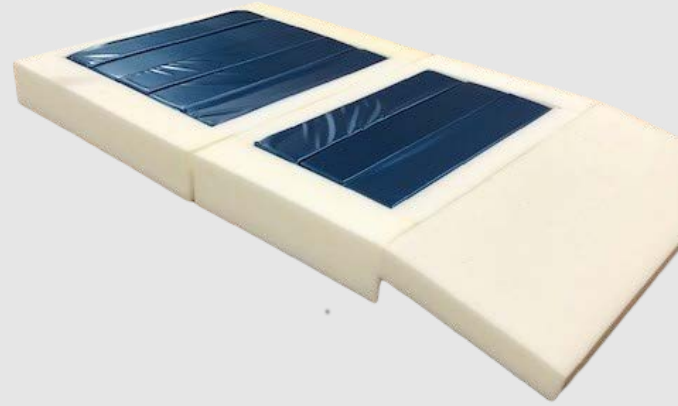

### Available on:




#### Deluxe

Premium Aesthetics with Powered Height with OneButtonSafeTurn, Night Light, and Battery BackUp



#### Standard

Powered Height with OneButtonSafeTurn

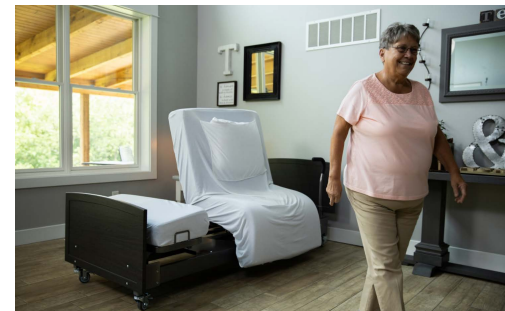

#### Fixed Height

Fixed Height with OneButtonSafeTurn

PROUDLY MADE IN THE U.S.A.

(877) 867-7365

**MED-MIZER.COM**



is at the forefront of innovating medical technology for post acute care. We are focused on partnering with you to help address healthcare challenges in the following key areas—




**FALLS PREVENTION**




**BARIATRIC CARE**



**MOBILITY**

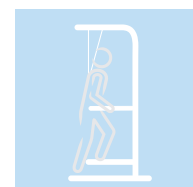

**REHABILITATION**

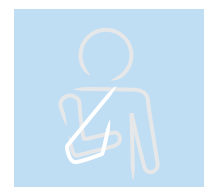

**WOUND CARE**